



ramcoTM
value of experience

Health & Home Care



GYNECOLOGY

Surgical Drapes And Pack



DESIRE TO CHANGE

From Reusable To Disposable




*Disposable drapes avoid repeated steps.
Save resource and offers safety.*

Washing & Cleaning



Drying, Folding and Packing

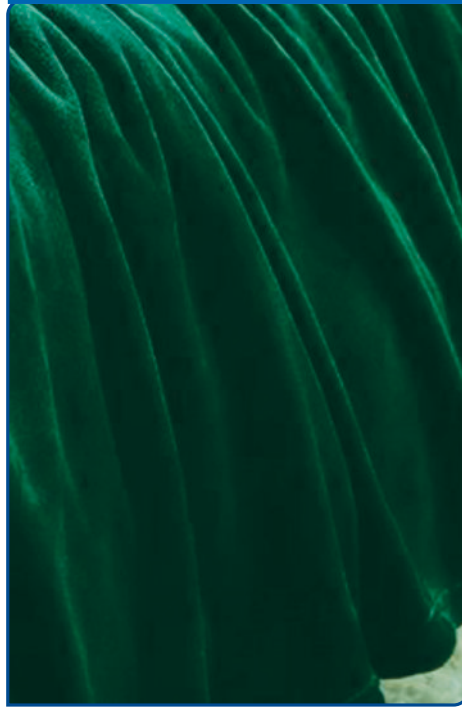

Sterilization



LINT is Source of Infection And Allergy



Blood Stains, Chemical Stain & tears may lead to discomfort



ADVANTAGES OF DISPOSABLE DRAPES OVER REUSABLE LINEN

Disposable Drapes avoid repeated Washing, Cleaning, Drying, Folding, Packing & Sterilization Step.
Disposable Drapes saves resources Water, Power, Time & Manpower
Disposable Drapes offers High Barrier to Blood and Body Fluids.
Disposable Drapes are cost effective.
Disposable Drapes are Lint free, Tear free - Hence Safer and Comfortable.

A COMPARATIVE STUDY

COMPONENT	SAFE WEAR	CHEAP ALTERNATIVE	SAFE WEAR ADVANTAGE
1. Incise film	High quality PU film	LD/PE films	<p>SAFEWEAR PU incise film is</p> <ul style="list-style-type: none"> ● Breathable : Allows air and moisture to transfer through the film keeping the covered area and comfortable. LD and PE film are non-breathable ● Conformable : Unlike LD/PE films, our incise film conforms easily to the shape of body contours ad creates a soft layer like a second skin allowing for cleaner and more precise incisions. Bacteria impenetrable : Creates a complete sterile field around the surgical site by creating a complete barrier to bacteria and other micro-organisms ● Low Glare : Unlike LD/PE films which create a glare through reflection of light, our incise film is extremely low glare in order to provide the surgeon with a clear, unobstructed view of the incise area during surgery
2. Adhesive	Skin friendly medical-grade adhesive	Industrial adhesive	<ul style="list-style-type: none"> ● The adhesive used by fine wear is carefully quality-controlled tested for cytotoxicity to ensure zero skin-irritation or skin allergies due to application ● The SAFEWEAR adhesives have the right balance of adhesive strength and peel strength to ensure proper adhesion during surgery and easy, pain-free peel-off after surgery ● Cheaper alternatives use normal industrial grade adhesive which do not care for the above properties and can cause skin problems
3. Fabric	High quality SMS / Sontara	Spunbond	<ul style="list-style-type: none"> ● SAFE WEAR uses only high quality SMS material, Which is a trilaminate fabric having 3 layers, while spunbond is a single layer cheap fabric normally used for making non-woven shopping bags. ● The middle layer of meltblown ensures over 98% bacterial filtration efficiency to minimize the chances of cross infection. spunbond, on the other hand, has almost zero bacterial filtration properties. ● Our SMS fabric has a high hydrostatic head to create a strong a liquidpenetration barrier, thus keeping the patient dry and comfortable. The hydrostatic head of sunbond is much lower by comparison. ● Our SMS fabric is also soft on the skin, unlike spunbond which has a much rougher feel.
4. Lining Film	Hygiene film	LD film	<ul style="list-style-type: none"> ● SAFEWEAR hygiene film is made of high quality virgin PE material, which gives it very good drapability and comfort compared to LD. ● Unlike some films that contain fillers, SAFEWEAR hygiene film is completely skin friendly. ● LD films make crackling sounds when moved or crumpled, which can be very distracting during a surgical procedure. Fine Wear hygiene films provide for extremely quiet handling, allowing the healthcare staff to concentrate other surgery without disturbance during the producer.






LSCS PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	3 nos.
2	Hand Towel	6 nos.
3	Cesarean Drape with Collection Pouch	1 no.
4	Under Pad	1 no.
5	Baby Receiving Sheet	1 no.
6	Baby Wiping Sheet	1 no.
7	Cord Clamp	1 no.
8	Wrap Sheet	1 no.
9	Cauty Bag	1 no.
10	Mayo Trolley Cover	1 no.
11	Trolley Sheet	1 no.

Benefits of Safe Wear LSCS Pack :

- Cesarean Drape comes with **wide collection pouch & two suction ports** - Avoids spill out & Offers sufficient fenestration area.
- Fenestration with **quality incise** - ensures easier access and minimal invasive
- Single Cesarean Drape **Covers patient & Screen Area.**
- Under pad comes with thick fluffy - offers higher absorbcency and allows fluid to pass quickly into.
- Under pad comes with Polypropylene back sheet protects against leakage.
- Cauty bag comes with 3 compartments.
- Standard labeling and arrow signs (→) on drape for unfolding - Ensures Faster and Comfortable draping

LSCS PACK

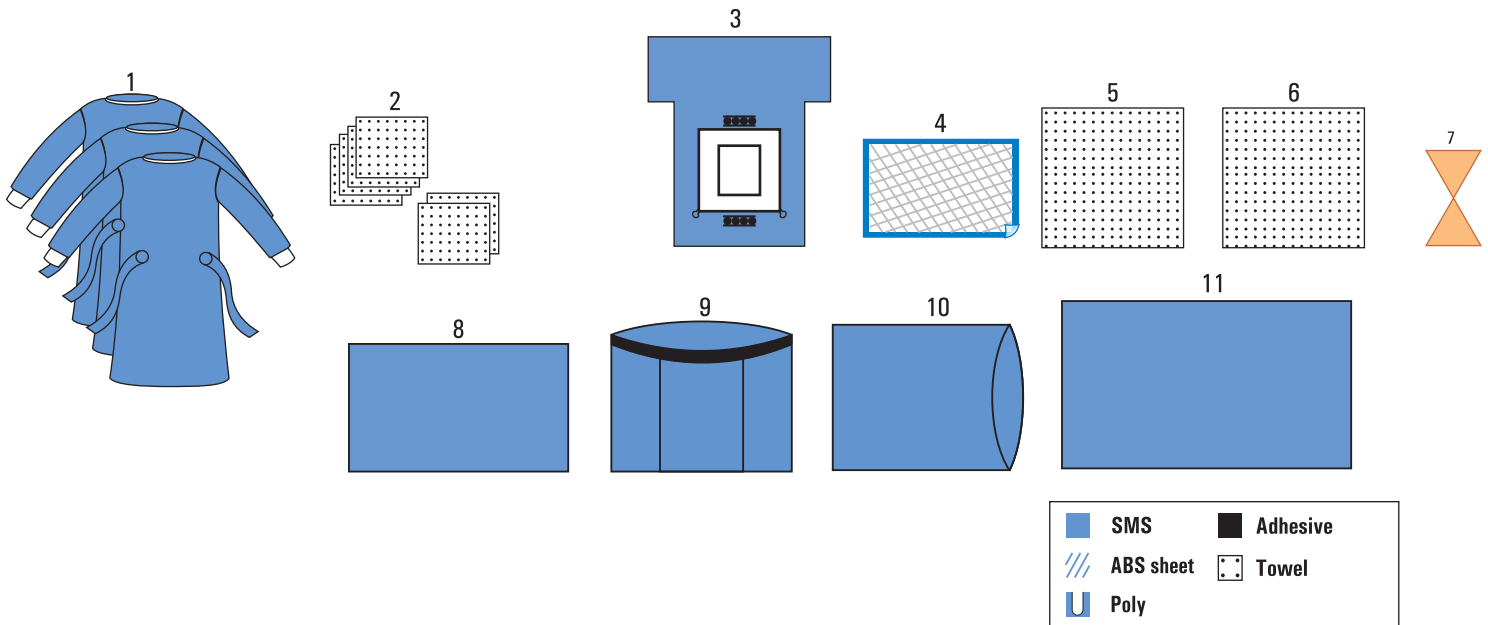


NORMAL DELIVERY PACK



Sl. No.	Description	Quantity
1	Surgeon Gown - Reinforced	2 nos.
2	Hand Towel	4 nos.
SAFE-1	Normal Delivery Drape	1 no.
4	Under Pad	1 no.
5	Baby Receiving Sheet	1 no.
6	Baby Wiping Sheet	1 no.
7	Cord Clamp	1 no.
8	Wrap Sheet	1 no.

SAFE-1

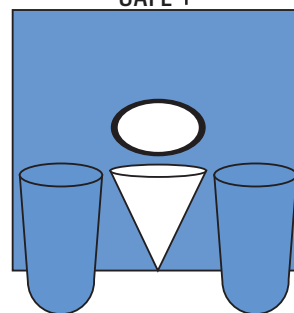

SMS	Adhesive Tape
ABS Sheet	Reinforced Sheet
Poly	


Benefits of Safe Wear Normal Delivery Pack :

- Normal Delivery drape comes with **wide collection pouch and suction ports - Avoids spill out.**
- Even Normal Delivery drape comes with leggings.
- Under pad comes **with thick fluffy - Offers higher absorbency and allows fluid to pass quickly into.**
- Under pad comes with Polypropylene back sheet protects against leakage.
- Standard labeling and arrow signs (→) on drape for unfolding - Ensures Faster and Comfortable draping.



ABDOMINAL HYSTERECTOMY PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	3 nos.
2	Hand Towel	6 nos.
SAFE-2	Large Drape with Adhesive	2 nos.
SAFE-3	Side Drape with Adhesive	2 nos.
8	Wrap Sheet	1 no.
9	Cauty Bag	1 no.
10	Mayo Trolley Cover	1 no.
11	Trolley Sheet	1 no.

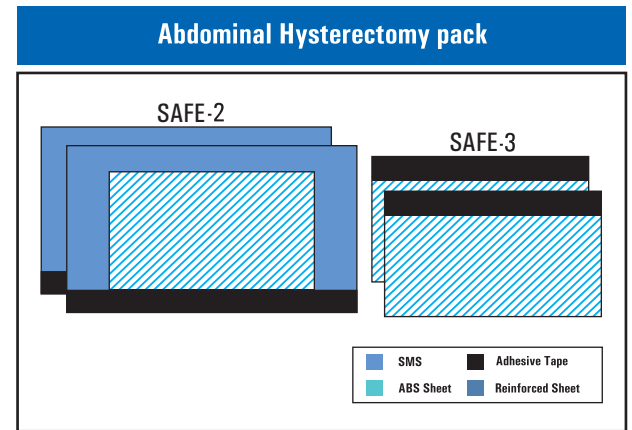


VAGINAL HYSTERECTOMY PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	3 nos.
2	Hand Towel	6 nos.
SAFE-1	Hysterectomy Drape	1 no.
SAFE-4	Large Drape with Out Adhesive	1 no.
4	Under Pad	1 no.
8	Wrap Sheet	1 no.
10	Mayo Trolley Cover	1 no.
11	Large Trolley Sheet	1 no.

Benefits of Safe Wear Vaginal Hysterectomy :

- Hysterectomy Drape comes with **wide collection pouch & two suction ports** - Avoids spill out and offers sufficient fenestration area.
- Even drape comes with leggings.
- **Under pad comes with thick fluffy** - offers higher absorbency and allows fluid to pass quickly into
- **Under pad comes with polypropylene back sheet** protects against leakage.
- **Standard labeling and arrow signs (→)** on drape for unfolding - Ensures **Faster and Comfortable** draping.

VAGINAL HYSTERECTOMY PACK

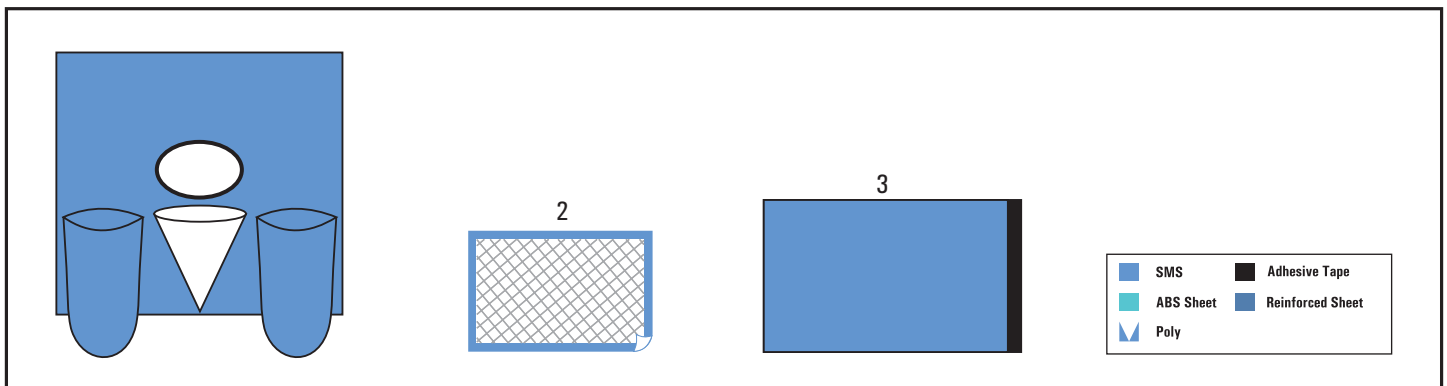


HYSTEROSCOPY PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	2 nos.
2	Hand Towel	4 nos.
SAFE-1	Hysteroscopy Drape	1 no.
8	Wrap Sheet	1 no.



IVF PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	2 nos.
2	Hand Towel	4 nos.
SAFE-5	Perineal Drape	1 no
SAFE-6	Medium Drape	1 no
SAFE-7	Leggings	2 nos.
8	Wrap Sheet	1 no.
12	Gauze Swab's	10 nos.

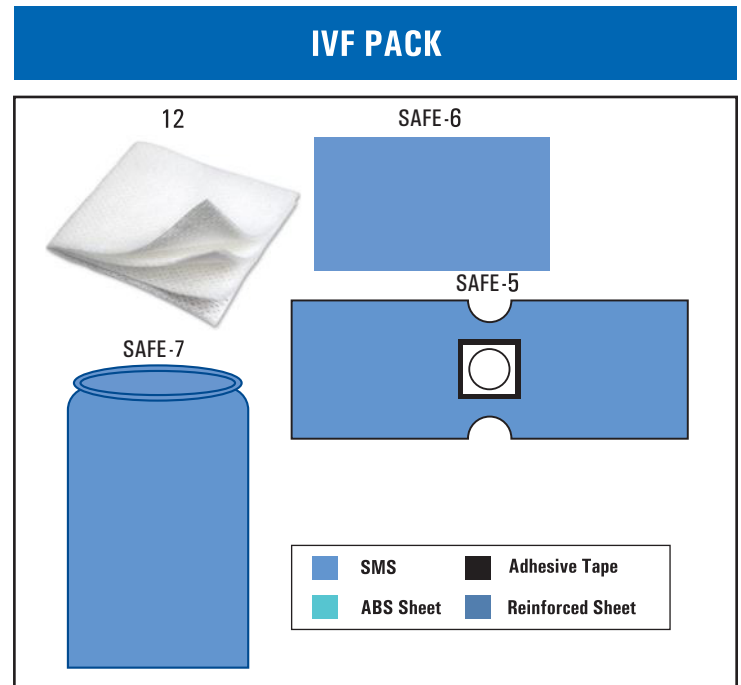

DILATION AND CURETTAGE (D&C) PACK

Sl. No.	Description	Quantity
1	Surgeon Gown	2 nos.
2	Hand Towel	4 nos.
SAFE-5	Perineal Drape	1 no.
SAFE-6	Medium Drape	1 no.
SAFE-7	Leggings	2 nos.
8	Wrap Sheet	1 no.

GYNECOLOGY - ACCESSORIES

Fine Wear Gynecology accessories are available in individual packs

Sl. No.	Description	Sl. No.	Description
1	Surgeon Gown	11	Trolley Sheet
2	Hand Towel	12	Gauze Swab's
3	Cesarean Drape w/ collection pouch	SAFE-1	Normal Delivery Drape
4	Underpad	SAFE-1	Hysteroscopy Drape
5	Baby receiving sheet	SAFE-2	Large Drape With Adhesive
6	Baby wiping Sheet	SAFE-3	Side Drape With Adhesive
7	Cord Clamp	SAFE-4	Large Drape with Out Adhesive
8	Wrap Sheet	SAFE-5	Perineal Drape
9	Cautery Bag	SAFE-6	Medium Drape
10	Mayo Trolley Cover	SAFE-7	Leggings
			Mopping Pad

DRESS FOR SURGICAL EXCELLENCE



CHANGE S DESIRE FOR EXCELLENCE

- ▶ Disposable drapes save resources and avoid time consuming processes for the hospital.
- ▶ Our medical fabric meets AAMI PB 70 (American) & EN 13795 (European) standards.
- ▶ Ensures faster & comfortable surgical draping.
- ▶ Committed to deliver quality products at affordable prices.
- ▶ Trusted product partner to more than 200 hospitals across india.
- ▶ Prompt delivery.
- ▶ Technical training to OT & CSSD personnel.

GET IN TOUCH

Regd. Office : RE-17/10, Raghunathpur Housing Estate
Raghunathpur, Kolkata - 700 059, WB, INDIA
Ph. - +91 9831 2322 47

@ ramcocal@gmail.com

www.ramcoindia.co.in

ISO : International Organization for Standardization (ISO), **AAMI** : Association for the Advancement of Medical Instrumentation, **ETO Sterilization** : Ethylene Oxide Sterilization, **CSSD** : Central Sterile Services Department, **SMS** : Spunbond-Melt-blown-Spunbond **LSCS** : Lower segment Cesaæan Section, **IVF** : In Vitro Fertilization, **D&C** : Dilation and Curettage.

Disclaimer : This is an independent publication owned by Fine Wear. Although great care has been taken in compiling and checking the information given in this publication to ensure that it is accurate, the author/s, the printer/s, the publisher/s and their servant/s or agent/s and purchaser/s shall not be responsible or in any way liable for any errors, omissions or in accuracies whether arising from negligence or otherwise howsoever, or for any consequences arising there from.



ramcoTM
value of experience
Health & Home Care